


New **Doublelo**

Product Manual

* WHEN OPERATING THIS PRODUCT, PLEASE READ
THE INSTRUMENT MANUAL CAREFULLY FIRST

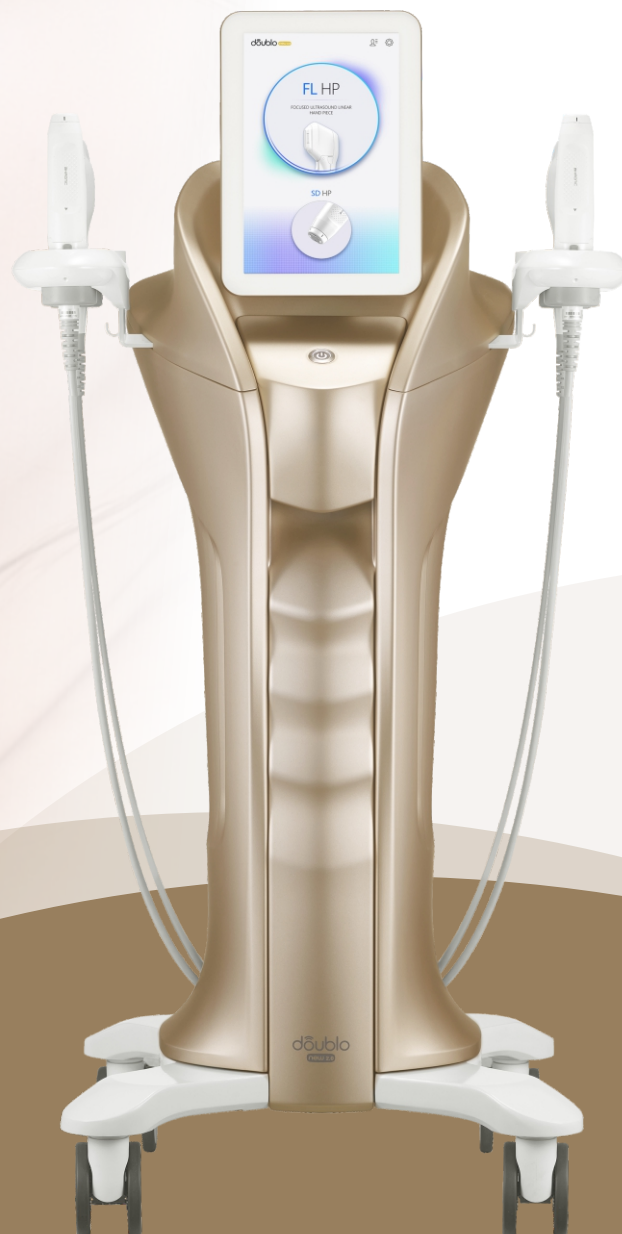

Preface

First of all, I would like to express my heartfelt thanks to the customers who purchased our instruments. This manual is prepared for customers to better use this instrument, in order to use the instrument correctly, to prevent damage to the operation and users, protect the instrument from damage, to avoid causing unnecessary property losses, please read the manual carefully before use. The Company considers itself responsible for the safety and reliability of the instrument only if the fitting, reassembly, improvement and maintenance are performed by personnel approved by the Company, the electrical equipment connected to the instrument conforms to national standards, and the instrument is used in accordance with the operating instructions. This is a kind of life beauty instrument. If the beauty organization using this instrument can not operate the instrument normally, it may cause the instrument function failure, and may endanger the health of the person. The Company will provide, for a fee upon request, a circuit diagram, along with a list of calibration methods and other information to assist the Customer in the maintenance of the parts of the instrument that the Company classifies as servable by the Customer. The product must be used in strict accordance with the usage method in this manual. The company shall not assume any responsibility for any consequences caused by misuse of methods other than the manual. The Company makes no warranties of any kind, including (but not limited to) any warranties of merchantability and fitness for a particular purpose. The Company shall not be liable for errors contained in this material or for incidental or consequential damages arising from the provision, actual performance or use of this manual. This user manual contains all the basic information of the instrument. This is the starting point for new users to start using this instrument, because it introduces the basic system, explains how it works, introduces the instrument installation and various operating methods, you can find the required information through the directory. This manual contains proprietary information protected by patent law. All rights reserved. Photographs, photocopies, reproductions or translations of any part of this manual into other languages are not allowed without the written consent of the Company. The contents contained in the manual are subject to change without notice to the user.

product list

1. Introduction to the instrument	
1.1 The principle of operation of the instrument.....	03
1.2 Appearance of the instrument.....	04
1.3 Interface Description.....	06
2. Notes on power on	
2.1 Handling safety.....	10
2.2 Electrical Safety	10
3. Operational Safety	10
4. Operating Instructions	10
5. Line drawing contraindications	11
6. Preparations before operation	13
7. Contraindications	13
8. Common faults and troubleshooting methods	13
9. Warranty Card	14

1. INSTRUMENT INTRODUCTION

1.1 HOW NEWDOUBLO WORKS

MMFU(Macro & Micro FocusedUltrasound) is "Macro & Micro Focused Ultrasound System" non-surgical facial lifting, firming and body contouring system. According to different skin types and treatment parts, the depth can be adjusted to treat the dermis and subcutaneous fat layers. The synergistic effect is produced by outputting energy at the same time with focused ultrasound + 4-pole radio frequency. The energy penetrates the inner and outer layers of the skin simultaneously without damaging the epidermis, resulting in a large amount of pain and redness.


1. SD dual technology synergy system of MFU (Micro Focus Linear Ultrasound) + RF (Radio Frequency)
2. Using specific ligament anchoring and pulling operation technology to achieve the effect of lifting, firmer, more adductor and more attractive skin.

1. Suspension lifting and tightening, painless and traumatic, effective at one time
2. Layered anti-aging septet, different depth probes
3. Three-dimensional improvement of skin age, skin texture and complexion
4. The original 2mm wrinkle-removing and tightening probe is aimed at the eyes and deep wrinkles - the MFU dual engine achieves a high peak value, increasing the energy by 3 times in one second
5. One lift and two firming, "ultrasonic D + heat M Ji" double guarantee, one machine can do it
6. -New Doublo ultrasound equipment can provide customized parameter setting services, a clear operation interface, and the status of the equipment is clear at a glance.

1.2 INSTRUMENT APPEARANCE AND INTERFACE DESCRIPTION





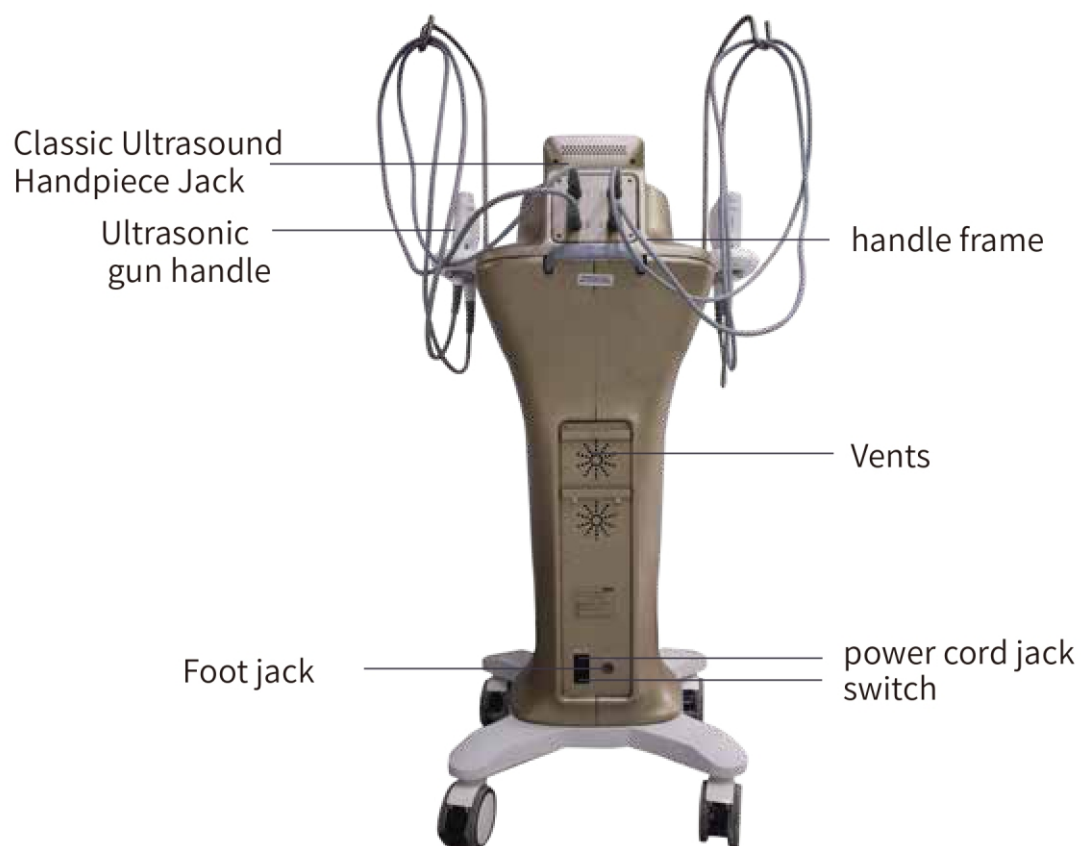



1.3 INTERFACE DESCRIPTION




Click to enter FL
Options interface

Click to enter SD
Options interface

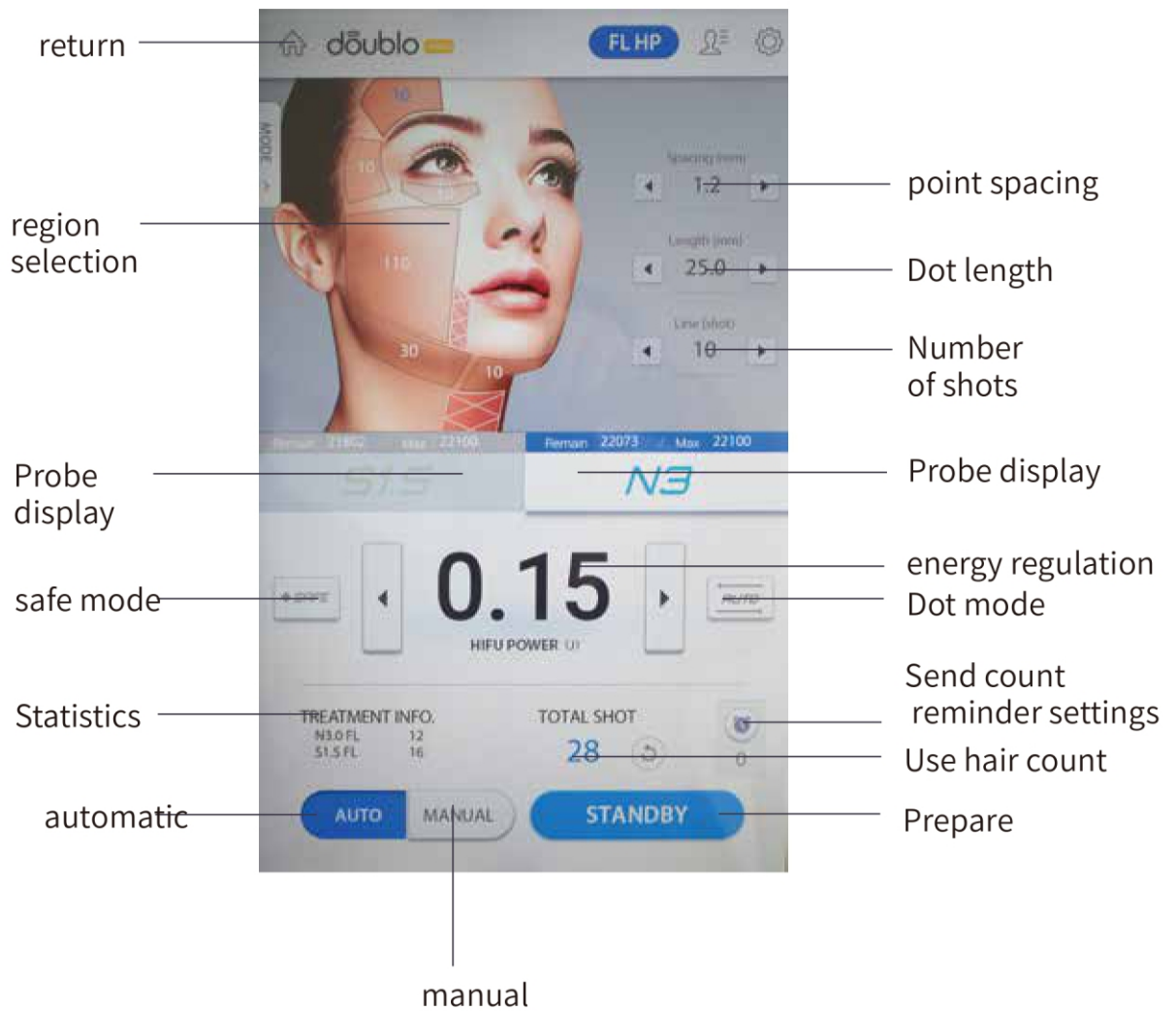

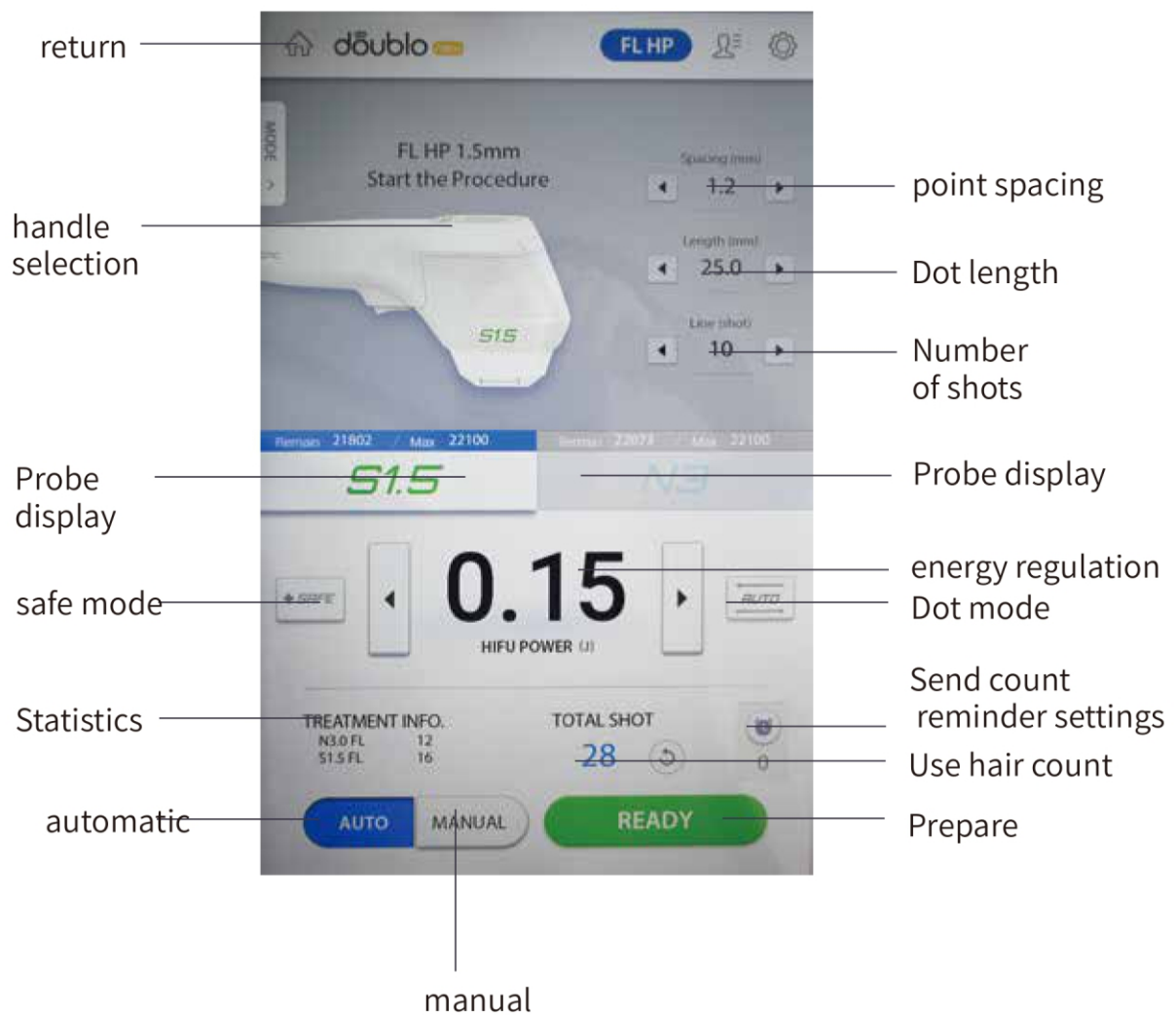

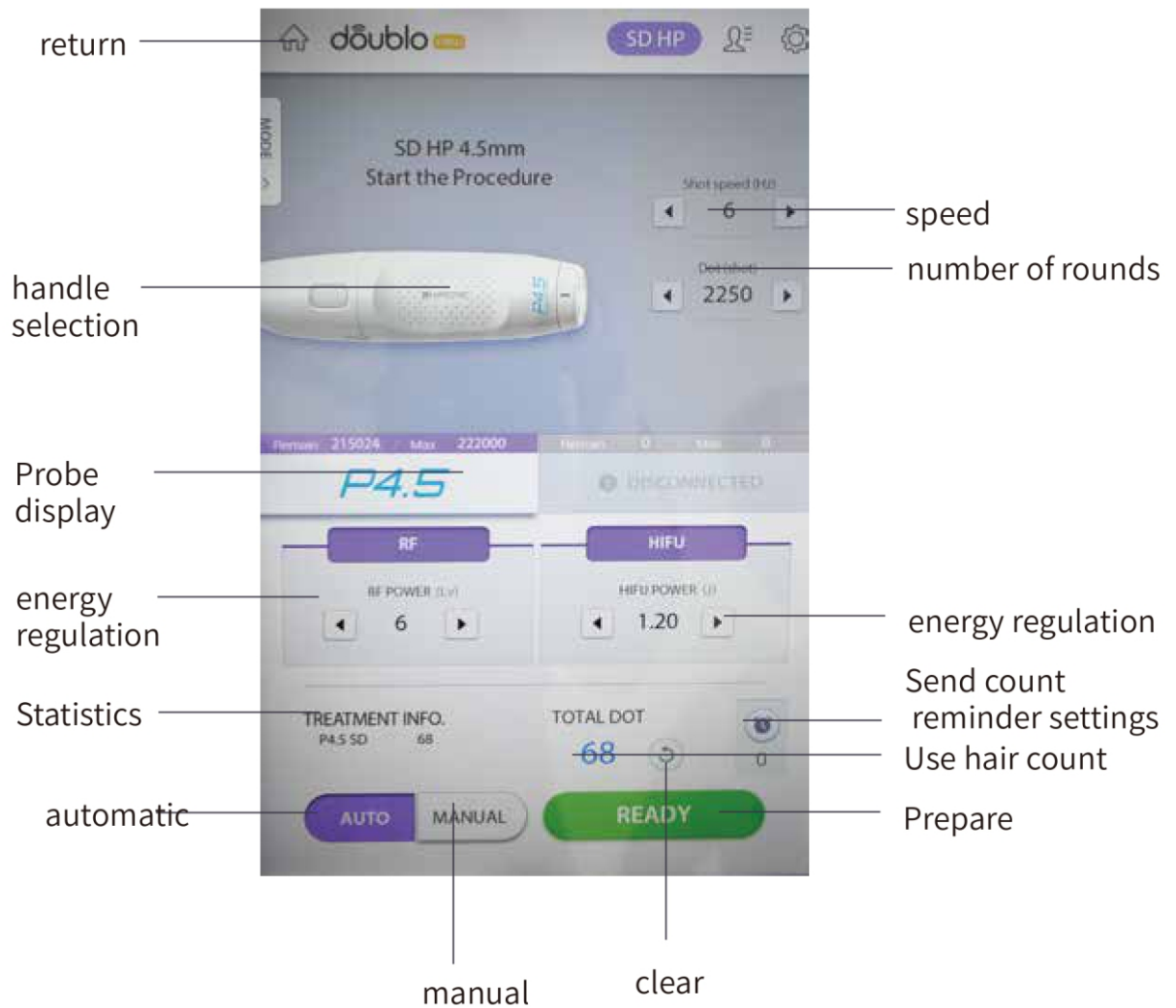

2. NOTES ON POWER ON

2.1 Safe handling

This operating system has precision instruments in memory, so users should try not to move this machine long distances. The shock-absorbing packaging provided by the original manufacturer must be used in the process, otherwise it may cause the optical system to be out of adjustment and affect the energy output.

2.2 Electrical safety

This operating system uses a single-phase 220V power supply, and the maximum output current is not more than 10A. In order to extend the service life of the machine and reduce the impact of power grid fluctuations on the machine, it is recommended that users refer to the use of a 3000w voltage stabilizer.

3. OPERATIONAL SAFETY

- 1) This instrument is dedicated for beauty salons and is not suitable for home users. The instrument needs to be operated by a professional beautician. Non-professional operation may cause unnecessary injury.
- 2) When you receive the instrument, please check all the accessories and check whether the machine is damaged. Refer to the list of instrument accessories. If the machine is damaged during transportation, please inform the factory.
- 3) The power supply used should be consistent with the specified power supply value marked on the machine, otherwise the machine may not work or even burn out the machine. This machine is grounded through the ground wire in the three-core power supply. Good grounding is essential for safe operation.
- 4) Please keep the original packaging of the instrument. When you need our assistance, you must wrap the machine in the original packaging before sending it for repair.
- 5) Do not disassemble the company's instrument or attempt to perform operations not described in this manual. All maintenance work must be performed by professional maintenance personnel approved by the company.
- 6) Do not install or operate the instrument near a water source or when your hands are wet, and be careful not to spill liquid on the instrument.
- 7) Do not place anything on the wire, and do not install the instrument in a place where the power cord is easy to be stepped on.
- 8) If the following situations occur, please unplug the instrument from the power socket and contact a professional service station.
 - a) Liquid splashed on the instrument, or the instrument was exposed to rain or water.


4. OPERATING INSTRUCTIONS

1. Do not operate on the bone area and nerve area;
2. Do not operate directly on the navel, posterior superior spine (PSIS), pelvis, ulna or sciatic nerve;
3. When caring for body parts, please treat on one side first, and then operate on the other side;
4. Avoid operating on specific parts that may cause tingling in the patient;
5. After stacking several times, pay attention to check the skin reaction and the sensitivity level of the patient. Continue treatment or reduce parameters according to the patient's tolerance level;
6. Ensure that the number of hairs is always evenly distributed over the marked treatment area;
7. The patient must select the appropriate probe according to the fat caliper measurement results of the intended operation site: FL-9.0/FL-6.0/SD-9.0

5. LINE DRAWING AND CONTRAINDICATIONS

TABOO AREAS AND FACIAL LINE DRAWING INSTRUCTIONS:

1. Mark the taboo areas: one finger beside the mouth, one finger beside the ear, the sunken temple, the corner of the eye, and one finger up on both sides of the brow;
2. Find the position of the jawline and draw a line along the edge of the jawline;
3. Find the highest point of the apple muscle, connect the nose wing to the apple muscle, and then connect the apple muscle to the middle of the ear;
3. Avoid 1-2cm at the corners of the mouth and connect vertically
4. Use the probe to operate the spacing, translate the spacing point, and connect the line;
5. When drawing the eye line, use your finger to touch out the position of the eye socket, and draw the line along the eye socket, and do not operate inside the eye socket. Hit the "X" position at the corner of the eye, do not operate;
6. The forehead starts from the taboo area of the brow, and moves outward in turn to draw lines.

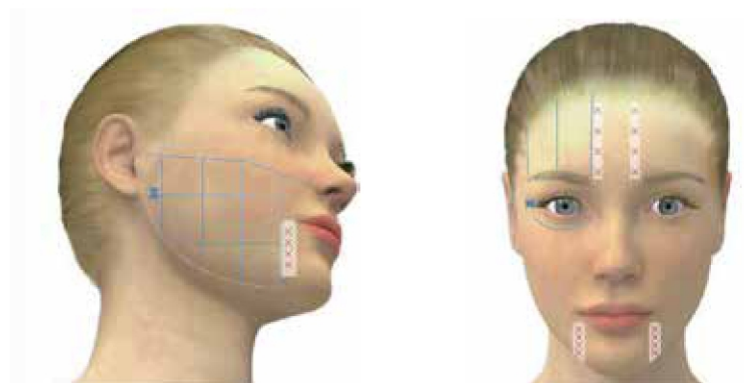

NECK LINE:

1. The mandible line is connected along the mandible, and then find the mandible near the throat and connect the line, as shown in the figure, the triangle area is connected vertically.
2. The neck first finds the position of the throat, marks the taboo area, and then translates and draws the lines in turn, and the opposite side is the same.



ABDOMINAL LINE:

1. Mark the taboo area first, open 2 fingers above and below the umbilicus, and draw a line along the lower edge of the rib. The ilium edge is drawn along the bone line.
2. Draw a line ruler on the abdomen, and connect the lines according to the position of the points, the opposite side is the same.



THIGH LINE DRAWING:

1. Find the nerve site and mark the taboo area.
2. Use the body to draw a line ruler, mark the position of the point and connect the line, the opposite side is the same.



6. PREPARATIONS BEFORE OPERATION

1. Ultrasonic gel, gel bowl and spatula
2. Cleansing basin, water, cleansing towel, cotton pad, cotton swab, paper towel
3. Disinfecting alcohol
4. Disposable towel and headgear
5. Disposable surgical gloves and masks

7. TABOO CROWD

The following groups of people are not suitable for nursing with this device:

1. Any wound on the facial skin
2. Having an electronic device or metal stand inside the body
3. Patients with wound and systemic disease
4. Patient has herpes virus infection
5. patients with autoimmune disease
6. patients with blood coagulation
7. pregnant or breastfeeding women
8. Patients with severe cystic acne
9. The wound has not been healed after surgery
10. Those receiving chemotherapy or cytoviral drugs
11. Patients with genetic allergies

8. COMMON FAULTS AND TROUBLESHOOTING METHODS

Failure phenomenon	Causes	Exclusions and solutions
Device not working	The power plug is not inserted into the socket The power switch is off Power fuse failure	Check power Check the power switch Replace the fuse
The handle does not work	Is the number of probes used up? The working head is damaged The interface on the back of the instrument is not installed correctly	Check the number of probes Replace working head re-install

9. WARRANTY CARD

1. The machine can be replaced with a new one within 7 days from the date of purchase (with the receipt) if there is a quality problem and the appearance is intact. One machine, the host is guaranteed for one year (except for wearing parts), other non-warranty scope or the warranty period has expired, also Responsible for repairs and updates, and appropriate repair fees - are required. If the instrument fails, it needs to be returned for repair. Within half a year from the date, the buyer shall bear the return freight, and the seller shall bear the freight; the freight for more than half a year shall be borne by the buyer.

2. Free warranty service will not be provided for the failures caused by the following users' personal reasons;

- ① Failure caused by unauthorized disassembly or modification of this product;
- ② Failure caused by careless use;
- ③ Failure caused by lack of reasonable maintenance;
- ④ Failures caused by failure to follow the correct instructions in the instruction manual;

Warranty Card

Name _____ phone _____ Postcode _____

Shipping address _____

product name _____ Purchase date _____

Maintenance records			
Return date	Fault content	Treatment plan	Maintenance company

Remarks: This card should be filled in by the buyer and seller. When returning the instrument for repair, this card must be filled out and sent with the instrument. After the instrument is repaired, the card will be sent back to the buyer along with the instrument. Please keep this card properly at that time.