



ULTRALIFT

When operating this product, please read the manual carefully first

Preface

First of all, we would like to express our heartfelt thanks to the customers who purchase our instruments. This manual is for customers to use the instrument better, To use the instrument correctly, Prevent damage to operations and users, To protect the instrument from damage, To avoid unnecessary loss of property, Please read the instructions carefully before using.

The Company is responsible for the safety and reliability of the instrument only if: The assembly operation, reloading, improvement and maintenance are carried out by the personnel approved by our company, The electrical equipment connected to this instrument conforms to the national standard, and the use of this instrument is carried out in accordance with the operation instructions.

This is a life beauty instrument.

If the beauty agency using this instrument can not operate the instrument formally, it may cause the failure of the instrument function and may endanger the personal health. The company will provide the circuit diagram at the request of the user, and also list the calibration method and other information to help the user to repair the equipment part of the company classified by the by qualified technicians.

This product must be used in strict accordance with the method of use in this manual, the misuse of methods other than the manual caused by all the consequences of their own responsibility, the company does not bear any responsibility.

The Company does not provide any form of security, including (but not limited to) real and suitable security against it for a particular purpose. The Company shall not be liable for any incidental or indirect damage caused by the error contained in this information or by the actual performance and use of this manual. This user manual contains all the basic information of this instrument. This is the starting point for new users to start using this instrument, because it gives a basic introduction to the system, explains how it works, introduces the installation of this instrument and various operating methods, can be through the directory to find the required information.

this manual contains proprietary information protected by patent law. All rights reserved. No part of this manual shall be photographed, copied, reproduced or translated into other languages without the written consent of the Company. The contents of the manual can be changed without notice to the user.

Directory

1. Introduction of Instrument.....	4
1.1 Instrument parameters.....	4
1.2 Appearance of instrument.....	5
1.3 Probe installation.....	5
2.Points for attention.....	6
2.1 Precautions before use.....	6
2.2 Contraindications.....	6
2.3 Attention in Customer Use.....	6
2.4 How to use the instrument.....	6
3.Introduction to the main interface of the instrument.....	7
4.Maintenance and repair.....	8
4.1 Maintenance.....	8
4.2 repair	8
5.Common Troubleshooting Methods.....	8
6.Manufacturer's liability.....	9
7.Storage and Transportation.....	9
8.Duration of use	9
9.List of accessories	10

1.1 Instrument parameters

Name	Parameters
input power	220V 50HZ
Output power	150W
Frequency	2-7MHz
Energy size	0.1-3.0
Machine colour	white

1.2 Instrument appearance





1.3 Probe installation

Instrument installation-install the handle head

- 1.Take out the handle head and open the protective film to ensure that the handle head interface is dry and free of foreign matter
- 2.Align the handle head interface with the handle, and the protrusions on the handle head correspond to the handle grooves
- 3 Gently insert the handle head horizontally into the handle, and the installation is complete when you hear(click)
- 4 When replacing or removing the handle head, press the safety lock on the side of the handle, and gently pull out the handle head horizontally .



2. Matters needing attention



2.1 Precautions before use

To prevent unauthorized use of the equipment, ensure that it is only used by trained and certified personnel.
In order to avoid possible injury during the nursing process, please do not operate near the eyes of the patient.
Avoid caring for the thyroid, thyroid cartilage and trachea.
Avoid direct nursing of breast tissue and main blood vessels.
If there is any failure or adverse event, immediately stop the nursing program.
Before taking care, remember to apply an appropriate amount of ultrasound gel in the care area.
Before outputting energy, check whether the probe is correctly and completely attached to the skin.

2.2 Contraindications (forbidden to be used in the following patients)

When pregnant or breastfeeding
Open wounds or injuries in the nursing area
Severe or cystic acne in the nursing area
Metal stents/electrical implants in the nursing area
Bioabsorbable mechanical implants
Any type of skin infection
Bleeding disease or dysfunction

2.3 Precautions for illness (Patients with the following diseases should pay attention to and confirm the appropriateness of care:

Anticoagulation therapy
Skin caused by autoimmune diseases, skin cancer and herpes simplex
Diabetes or epilepsy
Severe skin diseases
Exfoliation surgery or heart disease
Allergies or hyperplastic scars
Hemostatic dysfunction


System safety precautions

To avoid the risk of electric shock, be sure to check the probe, hand tools and power cord before use.
Do not use a power cord or probe that has been damaged or is leaking liquid.
Do not remove the cover or hand tool of the host, because it is not allowed to modify the device.
Do not lean on the device, and do not cause physical external shock to the device, as this may cause permanent system damage.
Do not use equipment with anesthetic gas or flammable materials in the surrounding area.

2.5 The steps of using the instrument

1. Before using the instrument, disinfect and wipe the probe used with alcohol, and wipe off the water droplets on the face;
2. Apply the product (photosensitive gel) evenly to the area to be treated;
3. Before the operation of the instrument, replace the required treatment head, and the energy range from low to high to the range that the experiencer can bear is better.
4. Keep the treatment head flat against the skin during the operation, and the operation spacing should be closely arranged;
5. During the operation, the experiencer should be constantly asked about the comfort level and the adjustment of the energy level to ensure that the experiencer accepts the experience in a pleasant state;
6. Focus on strengthening treatment in the problem area (increasing energy);

3. Introduction to the main interface of the instrument:

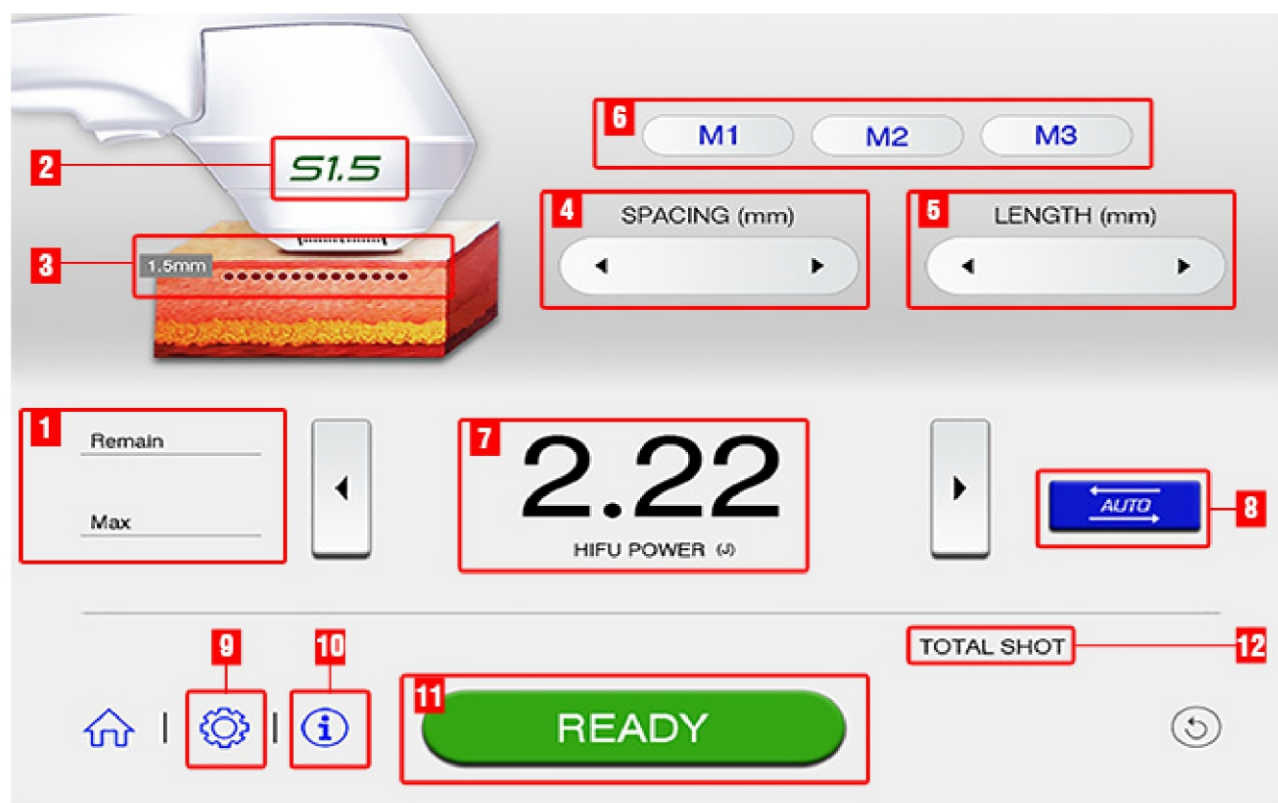

1. Remain remaining times, total times
2. Handle head information
3. Depth of the handle head, wavelength, spacing display (display only)
4. Set the distance between the condensation points
5. Set energy length (5mm-25mm, minimum adjustment interval 5mm)
6. M1, M2 save parameters (long press for 3 seconds to save)
7. Energy regulation
8. Repeat is continuous transmission mode Single is single transmission mode
9. Settings
10. Send number information
11. System status "standby" and "start"
12. the total number of times used

Step:

Preparation: Plug in the handle and probe, then adjust the dot (4) density and (5) intensity to start working (if you need to save the data, first adjust the required parameters, press and hold for 3 seconds (6) M1, M2 to save)

Start: first press (11) to start and then press the handle button, close to the skin, you can start treatment

Stop: first press the handle button and then press (11) to stop, switch the mode by pressing (6) to switch the mode.

Maintenance and repair

4.1 maintenance

4.1.1daily inspection

Management personnel must regularly inspect instruments and components.

Are power outlets reliably connected to the ground.

Please check whether the power cord is damaged before using the instrument, and stop using the instrument if it is damaged.

4.1.2 clean

Clean equipment regularly to ensure cleanliness and good operation of equipment.Wipe the device case with a wet cloth and use soapy water or a neutral cleaner if the device is powered off.

4.1.3replacement of fuse



When replacing the fuse, please turn off the power and pull the plug to prevent electric shock.

4.2 maintenance

Fter the purchase of equipment to fill in the warranty card, and send to the manufacturer or agent, to have the warranty effect;From the date of issuing the bill for the purchase of the machine belongs to the company's product quality problems mainframe free warranty for one year (note: with the purchase bill and warranty card warranty;This warranty is limited to the main machine, excluding spare parts and other consumables that need to be replaced periodically.After the warranty period, please contact the manufacturer directly to provide paid after-sales service.Maintenance of the machine must contact with the company's business department or dealers, do not repair with

5.Common faults and troubleshooting methods

The fault phenomenon	Cause of occurrence	Exclusion and resolution
Equipment not working	The power plug is not plugged into the socket The power switch is not on Power fuse failure	Check the power supply Check the power switch Replace fuse

Handle not working	The probe is not installed correctly The work head has been damaged The interface behind the instrument is not installed correctly	Check whether the probe is installed correctly Replace working head or silicone ring reinstall
--------------------	---	---

6. Manufacturer's liability

The products produced shall conform to the product standards. If the products do not conform to the product standards, the company is responsible for the replacement and return of goods.

The company shall bear the corresponding liability for personal and property damage caused by the manufacturing defects of the products, provided that the company shall not be liable for the following situations:

The purchaser or operation user fails to use the product in accordance with the operation specification in the instruction manual, or fails to use the product in accordance with the safety precautions and warning instructions.

Violation of regulations on transportation, installation, use, maintenance and storage of products, causing damage to products, human life or property.

The company carries out lifelong maintenance and technical service for the products, and will repair for users free of charge any manufacturing quality problems occurred during the stipulated quality assurance period.

When the accessories need to be replaced, the original accessories must be purchased from the company, otherwise the user shall bear the adverse consequences.

7. Storage and transportation management

Storage: the equipment should be stored in the original packaging box, in a well-ventilated room, the packing box should be high, the ambient temperature is -10°C~40°C, the relative humidity is not greater than 93%, not allowed to be placed in the vicinity of harmful gases, flammable, explosive substances and corrosive gases.

Transportation: the instrument under the packing condition is suitable for highway, railway, air and water transportation. In the process of loading, unloading and transportation, it shall prevent violent vibration and impact, and shall not be affected by moisture, and shall not be mixed with inflammable and corrosion-prone substances, mixed transport, specific requirements in accordance with the provisions of the order contract.

8. Use period


The warranty period of this instrument is 1 year (except for man-made damage).

The warranty period of spare parts is 3 months (except for man-made damage).

If the instrument beyond the service life continues to use, it shall be tested and approved by the engineer, and the performance of the instrument shall be tested regularly by the engineer in the subsequent use.

Scrapped instruments and accessories that have expired should be disposed of according to local regulations.

9.List of accessories

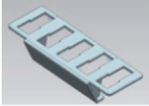Serial number	name	number	picture	Remarks
1	The host	1 Taiwan		
2	handle	2pcs		
3	Handle head	1.5、 2.0、 3.0、 4.5、 6.0、 9.0 1.5、 3.0、 4.5、 9.0		
4	The power cord	1article		
5	Film	14pes		
6	Acrylic sheet	2.5MM12 pes 3.5MM 6pes		
7	Handle bracket	2		
8	Probe holder	1		
9	Screw	10PES		
10	The instructions	1 this		