



UNILOG
DERMA SOLUTIONS



2In1

Diode Laser + Nd YAG Laser
Multifunctional Machine



Diode Laser

- Equipped with 5 removable treatment heads to meet the needs of different parts of the body
- Smart touch screen, easy to operate



1



Φ 8mm



5*5mm






10*10mm



12*20mm



20*40mm

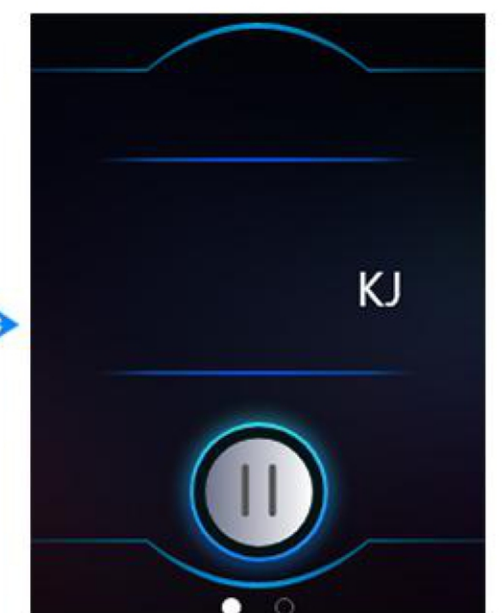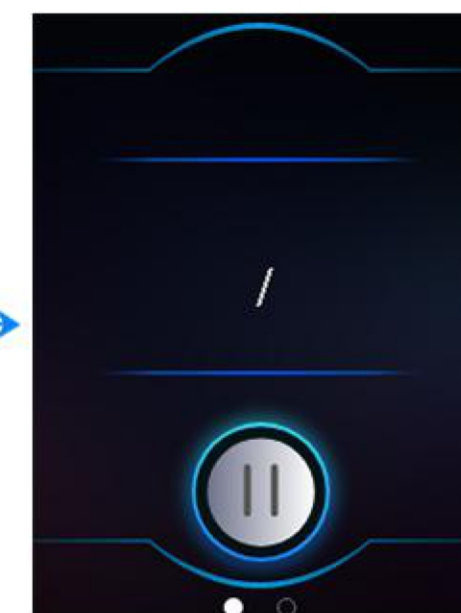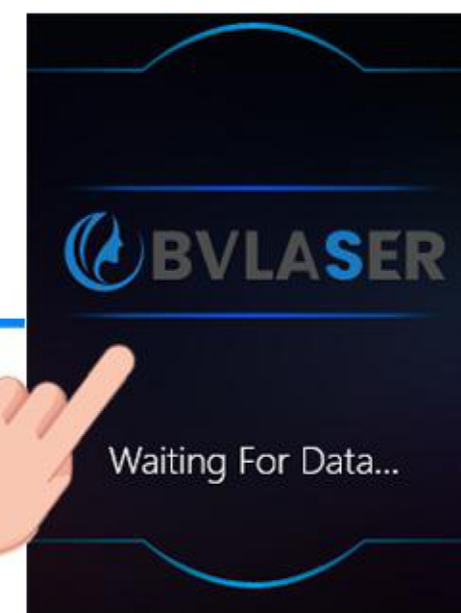



Larger spot size, Higher laser power

US Coherent imported laser device, high durability


Long service life



Spot size
20mm * 40mm



Four wavelengths 755nm, 808nm, 940nm and 1064nm laser, perfect for all hair types (thick, thin, black, gray, blond, brown, etc.) and all skin types (I, II, II, IV, V, VI)

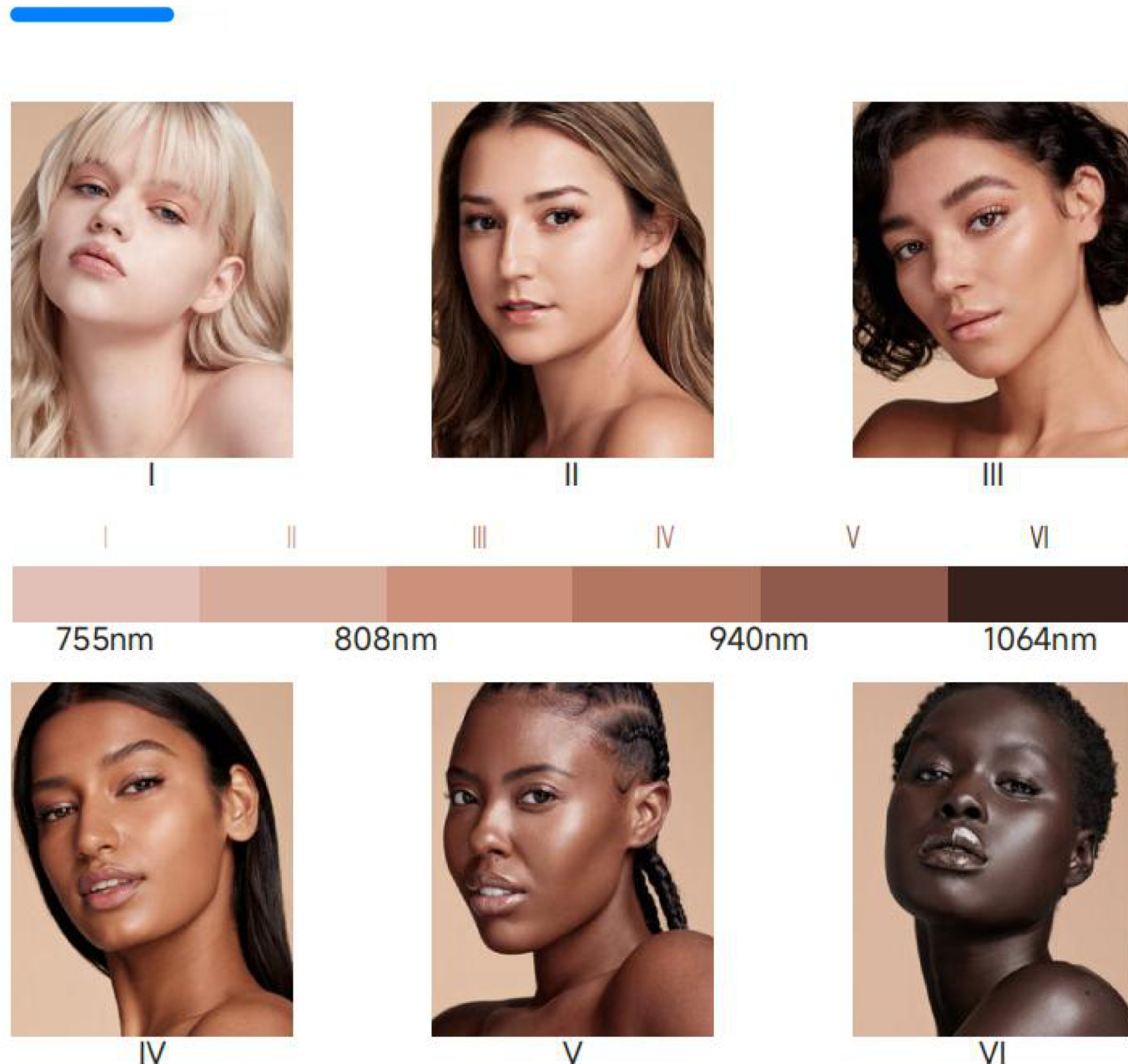

- 755nm** The Alexandrite 755nm laser is recommended for treating the fine and light colored hair.

With more superficial penetration, it targets the bulge of the hair follicle and is especially effective for the superficially embedded hair in some areas such as the eyebrows and upper lip.

- 808nm** The 808nm Diode laser has a deep penetration capacity that allows to treat the widest variety of hair types.

It focuses on the bulge and bulge of the hair follicle while moderate tissue depth penetration makes it ideal for treating the arms, legs, cheeks and beard.

- 940nm** 940nm is a useful wavelength for the treatment of coarser thicker hair on the skin.

It is effective in penetrating deeper into the skin, reaching deeper layers where hair follicles and subcutaneous adipose tissue are located.

- 1064nm** The Nd:YAG laser has a wavelength of 1064nm, focusing on treating darker skin. It can penetrate the deeper hair in the skin and has a greater water absorption which generates a higher temperature, increasing the temperature within the hair follicle for a more effective hair removal.

Q-switched Nd YAG Laser

Q-switched Nd YAG laser is available in different wavelengths for all skin types and can remove tattoos and pigments of all colors.



532nm

Remove light-colored pigments, such as shallow red, brown, pink, etc.



755nm

Remove melanin, improve skin tone, rejuvenate skin, and whiten skin, etc.



1064nm

Remove deep-colored pigments, such as black, blue, scarlet red, deep coffee, etc.



1320nm

Professional for skin rejuvenation, pimples removal, skin whitening, etc.



Before and After



Outstanding Features

THE MOST FULLY FUNCTIONAL LASER PLATFORM

The world's first 2-in-1 laser beauty platform integrates multiple functions, meeting almost all requirements of a beauty clinic.

SUPER STABLE

Use quality accessories and materials. Independent circuits and water paths ensure the stability of the laser platform.

COST-EFFECTIVE

Compared to two single-function machines, this 2 in 1 laser platform saves a lot of space and shipping costs. It is the most cost-effective multi-function laser beauty platform on the market.

LASER WATER FILTER

Double filtration ensures high purity of water. The intelligent self-check system will remind you to replace the filter cartridge in a timely manner. No frequent deionized water change.

SIMPLE BUT PROFESSIONAL USER INTERFACE

15.6-inch large color touch screen. Simple operation. Multi-language customization support. 100,000 customer usage records memory. System self-test for easy maintenance. The rental model makes it easier to manage customers.

Specification

Power	3000W
Screen	15.6-inch color touch LCD display
Cooling	Water + fans + semiconductor
Voltage	220V/50HZ, 110V/60HZ
Wavelength	Diode laser: 808nm/755nm/940nm/1064nm
	Nd YAG: 532nm/755nm/1064nm/1320nm
Spot size	Diode laser: 20*40mm, 12*20mm, 10*10mm, 5*5mm, ϕ8mm
	Nd YAG: 2-8mm
Energy	Diode laser: 1~120J/cm²
	Nd YAG: 1mj~1600mj
Frequency	Diode laser: 1Hz~10Hz
	Nd:YAG: 1Hz~10Hz

