



"Rich PTP Dual Pulse Delivery System Is Designed To Maximize Photoacoustic Impact While Mitigating Cumulative Thermal Load On Surrounding Tissue"

LUMASPOT

**HIGH POWERED,
HIGH DEFINATION RESULTS**

**NEW STANDARD IN
(Q-SWITCH ND:YAG LASER)**

LONG PULSE MODE:


Utilizes Millisecond Pulse Durations To Facilitate Extended Thermal Diffusion, Allowing The Target Chromophores To Absorb Energy Over A Timeframe That Exceeds Their Thermal Relaxation Time, Thereby Inducing Controlled Photothermal Coagulation.

Hair Removal

Facial Rejuvenation

Vascular Treatment

Onychomycosis (nail Fungus)

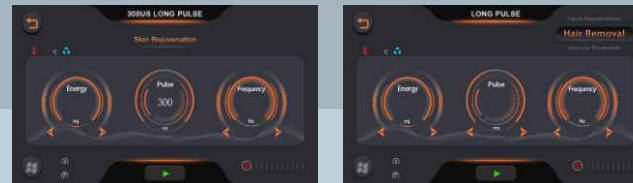

TECHNICAL SPECIFICATION

LASER TYPE	ACTIVE Q SWITCH ND:YAG LASER
SCREEN	15.6 INCH COLOR LCD SCREEN
LASER OUTPUT WAVE LENGTH	532/1064NM, 585NM, 650NM, 755NM, 1320NM
HANDPIECES: 7 HANDPIECE	(532NM/1064NM, 755NM, 585NM, 650NM, 1320NM, LONG PULSE 15MM, FRACTIONAL HANDPIECE)
MAXIMUM SINGLE PULSE ENERGY	1200MJ
MAXIMUM PTP ENERGY	2400MJ
MAXIMUM MULTIPLE PULSE MODE ENERGY	3600MJ
LASER PULSE WIDTH	5-10 ns
LASER RECURRENCE RATE	1-10HZ
AMING LIGHT	RED INFRARED DIODE LASER LIGHT
SPOTS DIAMETER	2-10MM(CONTINUOUSLY ADJUSTABLE)
CASE MATERIAL	ABS
LASER POWER SOURCE	220V/5A OR 110V/10A

LUMASPOT

**HIGH POWERED,
HIGH DEFINATION RESULTS**



📍 Sadhana Industries, A1, Swami Vivekananda Rd,
Opp. Tata Motors, Oshiwara, Jogeshwari West,
Mumbai, Maharashtra 400102

✉ info@unilogderma.com

🌐 www.unilogcaredevice.com





“Redefining The Spectrum of Skin Revitalization”

Is Next Generation Photomechanical Innovation With Advanced Modes Of Flawless Skin

LUMASpot


Rich Picosecond Technology:

Utilizes ultra-short pulse durations to induce Photoacoustic Shockwaves, shattering targeted chromophores into microscopic particles through mechanical stress rather than thermal accumulation.



Advanced PTP mode:


PTP Mode Optimizes Pigment Fragmentation Via Dual-pulse Energy Delivery, Maximizing Clinical Efficacy While Minimizing Peripheral Thermal Damage.

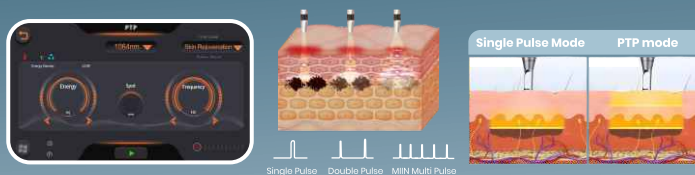

- Refractory Melasma (Dermal or Mixed)
- Dermal pigmentation on darker skin (Fitzpatrick IV–VI)
- Periorbital Pigmentation (Dark Circles)

Rich Multi Pulse Mode:

Quasi-Long Pulse or genesis mode Distributes high-density energy into a series of sub-microsecond pulses to induce **Selective Photothermolysis**, creating a controlled bulk heating effect in the dermis while preserving the epidermal threshold.

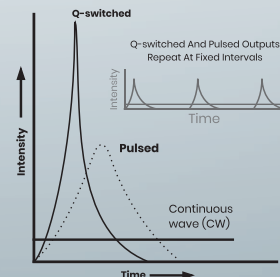

Carbon Peeling (Hollywood Peel)


BENEFITS :

Higher spot size, lower the cost to physician

- Safely treatment to all skin types
- Low maintenance ,minimum consumables required

User Friendly System Interface

- Multiple wavelength, multiple treatment modes.

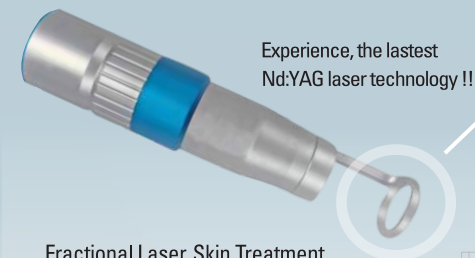


Experience, the latest Nd:YAG laser technology !!


Fractional Handpiece

Fractional Laser Skin Treatment


- ◆ High treatment efficacy by deep laser penetrating.
- ◆ Epidermis sparing during High fluence treatment.
- ◆ Minimize PIH, Hypo-pigmentation and side effect.
- ◆ 81 dots in 5 by 5mm spot size.

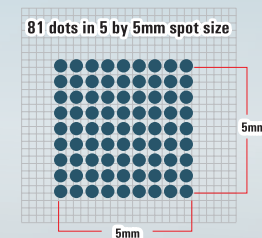

Top Hat Mode

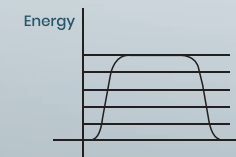

Beam Shape

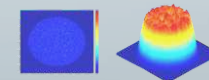

◆ Rejuvenation of Photoaged Skin

◆ Acne Scarring & Pitted Scars

◆ Striae (Stretch Marks)